

# Phillips Scientific

## Timing Discriminator

## NIM MODEL 715

### FEATURES

- UNEXCELLED TIMING CHARACTERISTICS
- FIVE TOTALLY INDEPENDENT CHANNELS
- 100 MHz OPERATION
- BOTH FAST VETO AND BIN GATE INHIBITING



### DESCRIPTION

The Model 715 is a high performance, five-channel constant fraction timing discriminator packaged in a single width NIM module. It exhibits typical time walk plus slewing of  $\pm 75$  pSEC from threshold to 100 times threshold.

Each channel has independent threshold adjustment, output width adjustment, and shaping delay. The threshold is variable from -25 mV to -1 Volt and is monitored from a front panel test point providing a DC voltage ten times the actual threshold setting. The outputs are non-updating with the output width adjustable from 5 nSEC to 150 nSEC. The shaping delay circuit is made complete by connecting the appropriate cable length between two LEMO connectors, allowing the user to optimize time resolution by easily matching the characteristics of each detector. A monitor is also provided to observe the constant fraction shaped pulse to verify that the delay is correct.

A fast veto input allows simultaneous inhibiting of all five channels to reject unwanted events early in the system. Similarly, a bin gate will inhibit the entire module when applied via the rear connector.

The 715 has three individually driven, current switching outputs per channel. Two are normal NIM levels and one complemented. They have typical risetimes and falltimes of 1 nSEC helping to preserve the excellent timing characteristics. The output transition times and shapes are unaffected by the loading conditions of the other outputs.

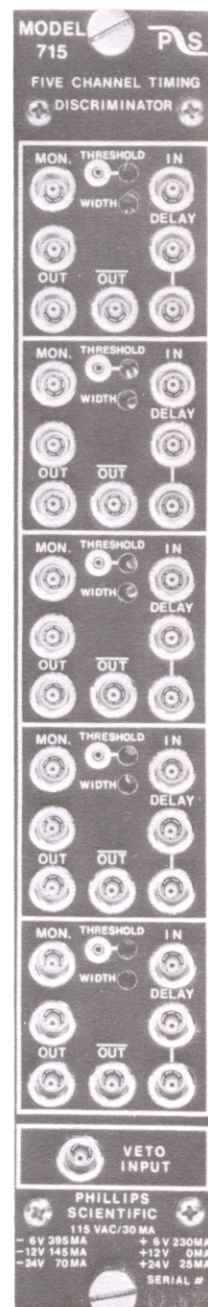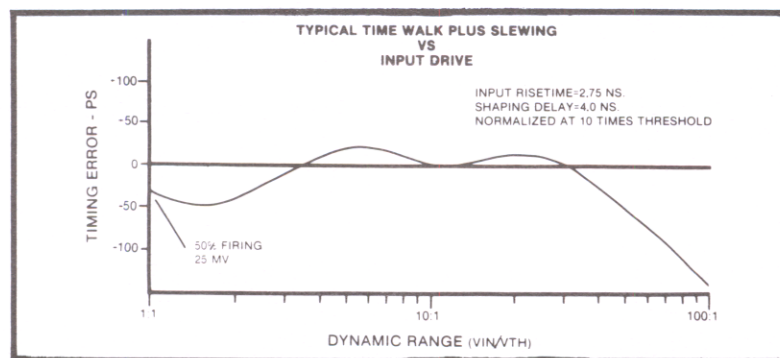

0983

**Phillips Scientific**

"A THEORY DEVELOPMENT COMPANY"  
150 Hilltop Road • Ramsey, NJ 07446 • (201) 934-8015 • Fax (201) 934-8269



## INPUT CHARACTERISTICS

### General:

One LEMO connector per channel; 50 ohms  $\pm 2\%$ , direct coupled; less than 5% input reflection for a 2.0 nSEC input risetime; input protection clamps at +7 V and -6 V and can withstand  $\pm 2$  amps ( $\pm 100$  Volts) for a duration of 1 uSEC with no damage to the input.

### Threshold:

Continuously variable from -25 mV to -1 Volt, 15-turn screwdriver adjustment; better than 0.2%/°C stability; front panel test point provides a DC voltage ten (10) times the actual threshold setting.

### Fast Veto:

One LEMO connector input common to all five channels; accepts normal NIM level pulse (-500 mV), 50 ohms direct coupled; must overlap the negative going edge of the input pulse plus the shaping delay time to inhibit all channels; 5 nSEC minimum input width.

### Bin Gate:

Rear panel slide switch enables or disables bin slow gate in accordance with TID-20893. Inhibits entire module within 10 nSEC from application of bin gate.

## GENERAL PERFORMANCE

### Shaping Delay:

Requires a 50 ohm coaxial cable; recommended delay range of 500 pSEC to 100 nSEC; the maximum delay is limited only by the cable attenuation, a factor of two attenuation can be tolerated without significant degradation of the time resolution; stability is better than 10 pSEC/°C. The shaping delay time should approximately equal the input risetime plus one nanosecond.

### Continuous Repetition Rate:

Greater than 100 MHz; with output width set at minimum. (1 nSEC shaping delay)

### Pulse Pair Resolution:

Better than 10 nSEC, with output width set at minimum. (1 nSEC shaping delay)

### Input to Output Delay:

Less than 10 nSEC; with 1 nSEC shaping delay.

### Multiple Pulsing:

One and only one output pulse regardless of input pulse amplitude or duration.

### Power Supply Requirements:

-6 V @ 390 mA	+6 V @ 240 mA
-12 V @ 150 mA	+12 V @ 0 mA
-24 V @ 45 mA	+24 V @ 35 mA

115 VAC @ 30 mA

**NOTE: All currents are within NIM specification limits permitting a full powered bin to be operated without overloading.**

### Operating Temperature:

0°C to 70°C ambient.

### Packaging:

Standard single width NIM module in accordance with TID-20893 and Section ND-524.

### Quality Control:

Standard 36-hour cycled burn-in with switched power cycles.

## OUTPUT CHARACTERISTICS

### General:

Three (3) LEMO connector outputs per channel; two normal NIM level outputs and one complementary output. The normal outputs deliver pulses of -16 mA (-800 mV across 50 ohms). The complemented output is quiescently -16 mA (-800 mV) and goes to zero mA (0 Volts) during output. Output risetimes and falltimes are less than 1.5 nSEC from 10% to 90% levels.

### CF Monitor:

One LEMO connector output per channel; drives 50 ohm load; permits observation of the shaped constant fraction pulse to verify the shaping delay is optimized for the input pulse.

### Width Control

One control per channel; 15-turn screwdriver adjustment; output width is continuously variable from 5 nSEC to 150 nSEC; better than  $\pm 0.2\%/^{\circ}\text{C}$  stability; nonupdating output regeneration.